

Quick User Guide to Comfort Audio FM system - Body worn TA version

Components

Teachers
Microphone



Pupils receiver with
neckloop




STEP 1: Switch on teachers microphone

Press and hold the button until the green bar has fully illuminated on the screen.



A STEP 2: Fitting and switching on pupils receiver


Place the neckloop around hearing aid users neck and plug in to receiver as shown in diagram. Switch hearing aid to the T Position.

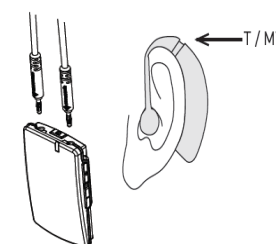

STEP 3 Wearing the microphone and receiver

Clip on microphone to clothing, the microphone works best around 30cm.

A necklace is provided if the clip is not possible to use.

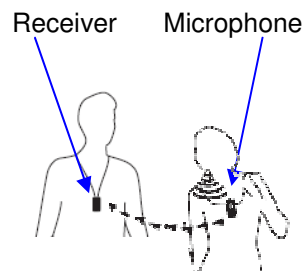


STEP 4 Increasing the volume

The volume of the neckloop can altered via the + and - switch on the side of the receiver.



STEP 5 Charging

Plug charger in to the mains using the charger provided.

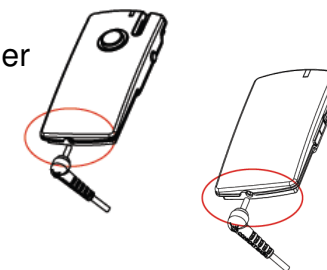

Connecting an audio source

Simply connect an audio jack from a laptop, pc, soundfield, etc in to the socket shown. This will transmit the sound to the hearing aid user.

